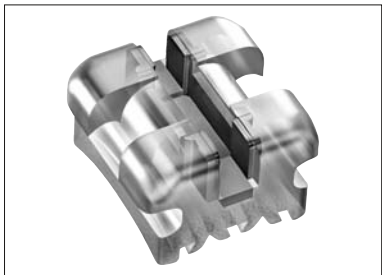
Elegance®

CE 0483

Gebrauchsanweisung
Instructions for use
Mode d'emploi
Modo de empleo
Modalità d'uso

D
DENTAURUM

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Instructions for use Elegance® brackets

Dear Customer

Thank you for having chosen a quality product from Dentaaurum. In order to use this product at its best for you and your patients, it is important to study and follow these directions for use carefully.

The written instructions cannot cover all eventualities during operation. For this reason our Hotline is available to answer any other questions and ideas that may arise.

Due to a constant development we recommend, even when you use the same products frequently, to study the enclosed latest issue of the instructions for use or refer to our website at www.dentaaurum.de.

Manufacturer

Dentaaurum GmbH & Co. KG · Turnstraße 31 · 75228 Ispringen · Germany

Working Steps

Inform the patient of the possible risks and adverse effects accompanying treatment with brackets, especially during debonding.

1. Tooth preparation

Check each tooth carefully for possible cracks or damage so that there are no complications during the debonding stage.

Please observe the manufacturer's instructions for the bracket adhesive you intend to use for conditioning the teeth.

2. Choice of bracket adhesive

Elegance® brackets can be bonded using any standard bracket adhesive. We recommend our own adhesives ConTec-, No-Mix-, Mix-adhesive and SmartBond®. Please note that it is absolutely necessary to adhere to the manufacturer's instructions for use of the relevant bracket adhesive, regarding the bonding of plastic brackets.

3. Removal and handling

Grasp the Elegance® brackets mesio-distally using bracket tweezers (eg. REF 025-277-00), remove from the packaging, and then observe the following working steps. The base of the brackets should not be touched by fingers or contaminated in any way, as this could have an adverse effect on the bonding quality.

4. Using the bracket adhesive

Prepare the bracket adhesive according to the relevant manufacturer's instructions. To ensure a secure and reliable bond of the adhesive on the plastic bracket base it is essential to apply a primer or conditioner to the base, which etches the base of the Elegance® brackets and additionally supplies a chemical bond between the plastic bracket and the adhesive. For the correct product please see the instructions for use of the adhesive manufacturer.

Use of light-curing adhesives, e.g. ConTec:

Apply a thin, even layer of ConTec primer for plastic brackets (REF 163-555-00 or 163-513-00) and polymerize for 20 sec. The ConTec adhesive must be applied immediately on to the bracket base.

Use of Mix- or No-Mix bonding systems:

Apply a thin, even layer of conditioner for plastic brackets (REF 163-500-00). The adhesive must be applied immediately on to the bracket base. The maximum drying time of 60 sec. should not be exceeded.

Use of Cyanacrylates, e.g. Smart Bond®:

No preconditioning of the bracket base is necessary.

Apply a small but sufficient amount of adhesive for covering the complete bracket base to the center of the base. This will guarantee an optimal bond.

5. Application and alignment of the brackets

After having applied the adhesive, place the Elegance® bracket onto the surface of the tooth in the customary manner.

It is essential to avoid adverse contact due to the risk of abrasion and/or premature bracket loss. If necessary a bite plate can be used in order to raise the bite. Remove all excess adhesive material from the bracket base before it sets. The removal of excess adhesive material is especially important to facilitate debonding after the treatment has ended.

6. Securing the arch wire

After the adhesive has polymerized (please see manufacturer's instructions for use) ligate as usual. For the initial treatment we recommend rematitan® „LITE“ or Tensic® arch wires.

7. Bracket debonding

The good bonding characteristics of Elegance® brackets means the process of debonding requires special attention: When debonding Elegance® brackets it is advisable to use a Weingart plier (eg. REF 003-120-00) or in the molar region, an angled bracket removing plier (REF 004-347-00). Grasp the bracket mesio-distally on the tie wings, and twist toward the left or the right without pulling away from the tooth.

8. Composition

The composition is included in the materials list, refer to the catalog or www.dentaaurum.de.

9. Safety information

The product should not be used, if there is a known allergic reaction to one or more of the material components.

10. Information on single-use products

The product is intended for single use only. Reconditioning of a product that has already been used (recycling) and its reuse on a patient is not permitted.

11. Additional Tips

- Technical information such as torque, angulation, size etc. can be found on the relevant product label or in the catalog information.
- Store Elegance® in sealed packaging.
- There are no particular precautions to be taken when disposing of Elegance® brackets after treatment.

12. With reference to our quality standards

Dentaaurum ensures a faultless quality of the products manufactured by us. These recommendations are based upon our own experiences. The user himself is responsible for the processing of the products. Responsibility for failures cannot be taken, as we have no influence on the processing on site.

13 . Explanation of symbols used on the label

REF	Order number
	Batch description
	CE Mark MP Class 2a
	Instructions for use available in the internet
	Not for reuse
	Codablock barcode
	Matrix code

Date of information: 02/10
Subject to modifications

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